



» Reusable surgical cutting instruments«





TEKNO-Medical Optik Chirurgie GmbH
Sattlerstr. 11
D-78532 Tuttlingen
GERMANY
SRN: DE-MF-000005822

Phone: +49 7461 17 01 0
Fax: +49 7461 17 01 50

Mail: mail@tekno-medical.com
Web: www.tekno-medical.com





Content

1	Scope	4
2	Inspections	4
3	Handling	4
4	Purpose	4
5	Indications	4
6	Contraindication	5
7	Combinations	5
8	Patient population	5
9	Disposal	5
10	Warnings	5
11	Instructions for reprocessing	5
11.1	On-site preparation	5
11.2	Transport	5
11.3	Preparation for decontamination	5
11.4	Manual pre-cleaning	5
11.5	Automated cleaning	6
11.6	Automated (thermal) disinfection	6
11.7	Functional testing, maintenance	6
11.8	Packaging	6
11.9	Sterilization	6
11.10	Storage	7
11.11	Information on the validation of the reprocessing	7
12	Additional Instructions	7
13	Reportable events	7
14	Warranty	7
15	Service and repair	7
16	Symbols	8
17	Product listing	8





In order to keep hazards to patients, users or, if necessary, third parties as low as possible, the instructions for use must be carefully observed. The application, reprocessing and testing of the instruments may only be carried out by trained specialists.



Reusable instruments from Tekno-Medical are delivered non-sterile and must undergo the complete reprocessing cycle (cleaning, disinfection and, if necessary, disinfection) before the first and each subsequent use. sterilization).

1 SCOPE



These instructions for use are valid for the reusable surgical cutting instruments of Tekno-Medical Optik-Chirurgie GmbH (see product listing in the last section).

2 INSPECTIONS

The instruments must be checked for functionality prior to every use.

Damages to the surfaces such as scratches, cracks, nicks, notches, etc., as well as bent parts mean that the instrument must not be used. The products must then be repaired by an authorized repair service or disposed of according to hospital standards.



Damaged products must not be used!

3 HANDLING

This instruction manual applies to instruments made of stainless steel, aluminium and titanium. Instruments made of pure titanium or titanium alloys should be treated in the same way as instruments made of stainless steel. There are no special instructions to follow.

Aluminium instruments must only be processed with non-alkaline, neutral disinfectants and cleaning agents and fully demineralized water, otherwise damage to the anodized surface can occur. Alkaline cleaning processes lead to colour fading and stains on colour-anodized surfaces after just a few cycles.

The instruments must not be subjected to excessive strain by twisting or levering, as this may result in damage or breakage of instrument parts. This product may only be used in medical facilities by trained healthcare professionals.

4 PURPOSE

Medical device with cutting or abrasive edges, with straight or curved blades for smoothing surfaces and/or cutting, separating or removing tissue such as skin, organs, bones or for removing foreign bodies., making incisions in the patient's body to allow the introduction of instruments into the surgical field.

The user decides on the basis of his expertise whether the instrument is suitable for the respective application.

Duration of use: Transient (intended for continuous use under normal conditions for a period of less than 60 minutes.) according to Regulation (EU) 2017/745

5 INDICATIONS

Instruments in this product group are used in almost all surgical procedures, e.g.

- Skin and tissue incision as well as other incisions with knives,
- making access to the surgical site,
- preparation of tissue and organ parts with scissors and clamps,
- Cutting and partial removal of bones or bone parts with saws, chisels, files, punches, bone spoons, also in the dental field,
- loosening of foreign bodies in or on the human body,
- Micro-scissors for fine, small procedures on tissue as well as on the eye,
- Snares / curettes: Separation of polyps of any kind with appropriate polyp lacers, snares, spoons, lancet needles for foreign body removal or sharp/blunt curettes for scraping out the uterus,
- Knives for amputation.





6 CONTRAINDICATION

In general, the use of the head support system is contra-indicated in cases where the use of other surgical techniques is indicated. In addition, there are contraindications:

- in case of unwillingness of the patient;
- if the technical requirements are not met.

Not for use on the central circulatory and nervous system within the meaning of the regulation.

The responsible physician or user must decide on the basis of the patient's general condition whether the intended application can be carried out.

7 COMBINATIONS

The products are not intended to be combined with or connected to other products.

8 PATIENT POPULATION

Apart from the contraindicated uses listed in these instructions for use, there are no restrictions on the patient population.

9 DISPOSAL

If the instruments can no longer be repaired and reprocessed, the instruments must be disposed of in accordance with the applicable country-specific regulations and laws.

Defective products must usually have gone through the entire reprocessing process before disposal.

10 WARNINGS

Instruments made of stainless steel must not be placed in physiological saline solution (NaCl), prolonged contact can lead to pitting or stress corrosion. Instruments may only be sterilized after prior cleaning and disinfection.

11 INSTRUCTIONS FOR REPROCESSING

In general, medical instruments may only be reprocessed by persons who have the necessary expertise for the intended activities. Detailed information on the preparation of instruments can be found in the "Red Brochure" of the AKI.

Under www.a-k-i.org you will also find links to laws, standards and reprocessing expert committees.

Due to the product design and the materials used, no defined limit of maximum feasible applications can be set. The lifespan of medical devices is determined by their function and careful handling. Frequent reprocessing has little impact on the product. End of product life is typically determined by wear and damage from use.

The readability of the labelling has been verified over 200 preparations.

11.1 On-site preparation

Immediately after use, remove coarse dirt from the instruments. Do not use any fixing agents or hot water (>40°C), as this will cause residues to freeze and may affect the success of cleaning.

11.2 Transport

Safe storage in a closed container and transport of the instruments to the reprocessing site to avoid damage to the instruments and contamination to the environment.

11.3 Preparation for decontamination

If possible, the instruments must be disassembled or opened for reprocessing.

The instruments must be stored on machine-compatible instrument carriers in a washer-safe manner. The condition of the instrument panels must not impair the subsequent cleaning and disinfection by means of ultra-sound or rinsing shadows.

11.4 Manual pre-cleaning

For manual cleaning / pre-cleaning, never use metal brushes, metal sponges or abrasive cleaning agents.

Place the instruments in cold water for at least 5 minutes. If possible, disassemble the instruments and clean them under cold water with a soft brush until no residue is visible. Pressure flush cavities, holes and threads with a water gun for at least 10 seconds (pulsed method, minimum pressure 2 bar).

Place instruments in an ultrasonic bath at 40°C for 15 minutes with 0.5% alkaline or enzymatic cleaner and sonicate. Remove instruments and rinse with cold water.

The cleaning solution should be changed at least once a day, more often if necessary. Too much contamination impairs the cleaning effect and increases the risk of corrosion.

National laws and guidelines must be observed.





11.5 Automated cleaning

Place the opened instruments in a sieve tray on the slide-in carriage and start the cleaning process. Disassemble instruments as much as possible.

Step	Parameter	
Pre-rinsing 1	Rinsing-temperature + water quality	Cold tap water
	Exposure time	60 s
Pre-rinsing 2	Rinsing-temperature + water quality	Cold tap water
	Exposure time	180 s
Cleaning	Cleaning-temperature	45 °C
	Water quality	Tap water
	Exposure time	300 s (worst case condition) / RKI Recommendation: 600 s
	Detergent	Neodisher Medizym
	Concentration	0,50 %
Neutralization	Rinsing-temperature	40 °C
	Water quality	Tap water
	Exposure time	180 s
	Neutralizing detergent	Neodisher Z
	Concentration	0,10 %
Post-rinsing	Rinsing-temperature	40 °C
	Water quality	Deionized water
	Exposure time	120 s

11.6 Automated (thermal) disinfection

Step	Parameter	
Thermal Disinfection	Disinfection-temperature	90 °C (A ₀ 3000)
	Water quality	Deionized water
	Exposure time	300 s
Drying	Drying of the outside of the instruments by the drying cycle of the washer-disinfector. If necessary, manual drying can also be achieved with the help of a lint-free cloth. Dry cavities and channels of the instruments with sterile compressed air. Allow products need to cool down to room temperature.	

For UK: The washer disinfection cycle should run at a minimum of 90 °C for a minimum of 1 minute.

11.7 Functional testing, maintenance

Visual inspection of cleanliness; If necessary, assembly of the instruments, maintenance and functional testing.

If necessary, repeat the reprocessing process until the instrument is visually clean.

Treat instruments with moving parts with care oil, e.g.: TK95 100-00. Only close instruments with locks in the first notch.

Defective or damaged instruments must be discarded immediately.

11.8 Packaging

Select standard-compliant packaging of the instruments for sterilization according to DIN EN ISO 11607-1, DIN EN 868-2 and DIN EN 868-8.

11.9 Sterilization

Sterilization of the products using a fractionated vacuum process (according to DIN EN ISO 17665-1), taking into account the respective national requirements.

Pre-vacuum	3 times
Sterilization temperature	134 °C
Sterilization time	5 min
Drying time	20 min.

The use of other sterilization methods is beyond our responsibility.

Sterilization parameters UK:

Sterilization of the instruments by applying a fractionated pre-vacuum process (according to Health Technical Memorandum 01-01 Part C):

- Heat up to a minimum sterilization temperature of 134° - 137°C (maximum temperature 137°C).
- Minimum holding time: at least 3 min.



11.10 Storage



The sterilized instruments must be stored in suitable packaging in a dry, clean and dust-free environment and at a constant level of humidity. The distance between the floor and the shelf should be at least 30cm. The storage period is to be determined by the user himself.

11.11 Information on the validation of the reprocessing

The following materials and machines were used in the processing process:

Detergent	Neodisher Medizym 0,5 % (v/v)
Neutralisator	Neodisher Z 0.1% (v/v)
Washer-disinfector (RDG)	Miele PG 8535
Steam-autoclave	Lautenschläger ZentraCert
For details see test reports: 23277 / 23278 / 23279 (CleanControlling Medical GmbH & Co. KG, 08-2021)	

12 ADDITIONAL INSTRUCTIONS



If the chemicals and machines described above are not available, it is up to the user to validate his process accordingly. It is the responsibility of the user to ensure that the remanufacturing process, including resources, materials and personnel, is suitable to achieve the required results.

The state of the art and national laws require the following of validated processes.

13 REPORTABLE EVENTS



In accordance with the requirements of the Medical Device Regulation EU MDR 2017/745 and our quality management system, even the smallest problems with this product should always be reported to TEKNO.

If you cannot reach us directly for reportable events, please send an e-mail to:

safety@tekno-medical.com

Serious incidents must also be reported to the competent authority in their place

14 WARRANTY

The products are made of high-quality materials and undergo quality control before delivery. However, if errors occur, please contact our service.

Tekno cannot guarantee that the products are suitable for the respective procedure. This must be determined by the user himself. Tekno accepts no liability for any incidental or consequential damages.

Tekno assumes no liability if it can be proven that these instructions for use have been violated.




Attention: In case of use of the instruments in patients with Creutzfeldt-Jakob disease or its variants (vCJD, BSE, TSE), Tekno-Medical declines any responsibility for reuse.

15 SERVICE AND REPAIR

Do not carry out repairs or modifications of the product on your own. Repairs only by authorized personnel.

Defective products must have gone through the entire remanufacturing process before being returned for repair. For returns, use our RMA request form and decontamination certificate.

Forms at: https://www.tekno-medical.com/de/service/repair_service

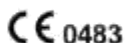

**16 SYMBOLS**

The symbols used in this instruction and on the label have the following meaning according to DIN EN ISO 15223-1:

	Attention!		Manufacturer
	Medical device		Date of manufacturing
	Non-sterile		Read the instructions for use
	Catalogue number		Protect from sunlight
	Lot-number		Store in a dry place
	Unique device identification		
	CE marking with number of the Notified Body mdc – medical device certification GmbH Kriegerstrasse 6, D – 70191 Stuttgart		

17 PRODUCT LISTING**REF**

Printed on: 19.03.2025

6500-08	8463-13-A	9718-17	24450-20	28003-05ZA	34470-25	40200-45	52675-13
6500-10	8463-13-B	9718-17KT*	24450-30	28003-06	34473-01	40202-10	52675-14
6500-12	8463-16	9718-17SR	24450-40	28003-06A	34473-02	40202-20	52675-15
6500-14	8463-16-A	9725-18	24451-01	28003-06Z	34475-01	40202-25	52676-01
6500-16	8463-16-D	9725-20	24451-02	28003-06ZA	34475-02	40202-30	52676-02
6500-18	8463-16-E	9725-23	24451-03	28004-01	34480-25	40202-35	52676-03
6500-28	8463-16SR	9729-28*	24451-04	28004-02	34485-20+	40204-01	52676-04
6500-30	8463-16SU	9731-23*	24451-10	28004-03	34485-25+	40204-02	52676-05
6500-32	8463-18	9822-23	24451-20	28004-04	34485-30	40204-03	52676-06
6500-34	8463-18SU	9827-22	24451-30	28004-05	34487-20	40210-01	52676-07
6500-36	8464-16*	9842-19	24453-04	28004-06	34489-19US	40210-02	52676-08
6500-38	8465-16*	9842-22	24453-05	28005-01	34489-24	40211-01	52676-09
6501-06	8469-12*	9847-21	24454-03	28005-01A	34490-00	40211-02	52678-01
6501-08	8477-24	9922-13	24455-00	28005-01ZA	34490-01	40212-01	52678-02
6501-10	8523-12	9922-13SU	24455-01	28005-02	34490-02	40212-02	52678-03
6501-12	8524-12	9923-13	24455-02	28005-02A	34490-03	40215-01	52678-04
6501-14	8525-12	9923-13SU	24456-00	28005-02ZA	34490-04	40220-01	52678-05
6501-16	8533-12	9923-13SUB	24456-01	28005-03	34490-05	40220-02	52678-06
6501-18	8534-12	9924-13 SU TC	24456-02	28005-03A	34490-06	40221-01	52678-08
6502-08	8535-12	9925-13 SU TC	24459-02*	28005-03ZA	34490-07	40221-02	52678-09
6502-10	8542-12	9936-08*	24459-20*	28005-04	34490-08	40225-01	52678-10
6502-12	8543-12	9942-14+	24460-01	28005-04A	34490-09	40430-01	52678-11
6502-14	8553-12	9982-10	24460-02	28005-04ZA	34491-00	40430-02	52678-12
6502-16	8562-11	9982-11	24460-03	28005-05	34491-01	40430-03	52678-13
6502-18	8562-12	9983-10	24460-04	28005-05A	34491-02	40455-00	52678-14
6505-08*	8563-11	9983-11	24460-05	28005-05ZA	34491-03	40460-00	52678-15





6506-10	8571-11	21101-15	24460-06	28005-06	34491-04	40475-00	52680-04
6506-12	8571-14	21101-19	24460-10	28005-06A	34491-05	40480-00	52680-05
7000-03*	8572-14	21101-21+	24460-20	28005-06ZA	34491-06	40760-09	52680-06
7000-10*	8573-14	21101-23	24460-30	28006-01	34491-07	40850-01	52680-07
7000-18*	8582-16	21102-15	24460-40*	28006-01A	34491-08	40850-02	52680-08
7005-01	8582-16D	21102-19	24461-01	28006-01Z	34491-09	40850-03	52680-09
7005-02	8587-16	21102-21+	24461-02	28006-01ZA	34491-20	40860-01	52680-10
7012-01	8589-16	21102-23	24461-03	28006-02	34491-20US	40860-02	52681-01
7012-02	8591-16	21103-15	24461-04	28006-02A	34491-25	40862-01	52681-02
7012-03	8592-12	21103-19	24461-05	28006-02-C	34491-25US	40862-02	52681-03
7012-04	8597-12	21103-21+	24461-10	28006-02Z	34492-01	40880-01	52681-04
7012-05	8602-16	21103-23	24461-20	28006-02ZA	34492-02	40880-02	52681-05
7012-06	8602-16SU	21104-15	24461-30	28006-03	34492-03	40880-03	52681-06
7012-07	8602-16SUB	21104-19	24461-40	28006-03A	34492-04	40880-04	52682-05
7012-08	8602-18	21104-21+	24462-01	28006-03-C	34492-05	40881-00	52682-06
7012-09	8602-18SU	21104-23	24462-02	28006-03Z	34492-06	40882-00	52682-07
7012-10	8603-14-A	21121-01	24462-03	28006-03ZA	34492-07	48120-00	52682-08
7100-04	8603-16	21121-02	24463-01	28006-04	34492-08	48120-01	52682-09
7101-04	8603-16-A	21121-03	24463-02	28006-04A	34492-09	48120-02	52682-10
7102-04	8603-16-E	21229-24	24464-01	28006-04-C	34492-10*	48120-03	52682-11
7105-08	8603-16SU	21232-21	24464-02	28006-04Z	34492-11*	48120-04	52682-12
7110-03	8603-16SUB	21233-21	24464-03	28006-04ZA	34492-20	48125-01	52684-05
7110-07	8603-18	21234-24	24464-04	28006-05	34492-20US	48125-02	52684-06
7110-71	8603-18SU	21235-24	24464-05	28006-05A	34492-25	48125-03	52684-07
7111-03	8603-18SUB	21236-01	24464-06	28006-05-C	34492-25US	48125-04	52684-08
7112-03	8603-20*	21236-02	24464-10	28006-05Z	34810-30	48125-05	52684-09
7113-03	8603-23SU	21236-03	24464-20	28006-05ZA	34810-35	48125-06	52684-10
7113-05	8612-12	21236-04	24464-30	28006-06	34811-30	48126-01	52684-11
7115-03	8614-12	21236-05	24465-01	28006-06A	34812-30	48147-00	52684-12
7116-03	8616-12	21236-06	24465-02	28006-06Z	34813-30	48147-01	52685-01
7125-13	8618-12	21236-07	24465-03	28007-01	34814-30	48147-02	52685-02
7500-13	8622-10	21236-08	24466-01*	28007-01A	35000-00	48147-03	52685-03
7500-16	8623-10	21236-09	24466-03*	28007-01Z	35000-01	48147-04	52685-04
7500-19	8632-11	24001-12	24466-04*	28007-01ZA	35000-02	48147-05	52685-05
7500-22	8633-11	24001-35	24468-02*	28007-02	35000-03	48147-06	52685-06
7528-05	8636-11	24002-03	24471-00	28007-02A	35000-04	48147-07	52685-07
7529-05	8636-11SU*	24002-05	24471-01	28007-02Z	35000-05	48150-02	52685-08
7550-12	8637-11	24002-07	24471-02	28007-02ZA	35000-06	48150-03	52685-09
7900-00	8637-11SR	24002-39	24472-00	28007-03	35000-07	48150-04	52685-10
7901-00	8639-11	24003-31	24472-01	28007-03A	35000-08	48150-05	52686-05
7902-00	8642-09	24004-23	24472-02	28007-03Z	35000-09	48150-06	52686-06
7904-00	8643-09	24005-27	24472-03	28007-03ZA	35000-10	48150-07	52686-07
7905-00	8652-09	24005-30	24472-04	28007-04	35000-11	48150-08	52686-08
7907-00	8652-10	24005-34	24472-05	28007-04A	35000-12	49921-15	52686-09





7907-01*	8652-10GL	24010-30	24472-06	28007-04Z	35000-20	50311-16	52686-10
7908-00	8652-10SU	24011-27	24473-00	28007-04ZA	35001-00	50313-16	52686-11
8023-10	8652-11	24014-28	24473-01	28007-05	35001-01	50315-17	52686-12
8023-11	8653-09	24015-24	24473-02	28007-05A	35001-02	50319-20	52687-01
8023-13	8653-10	24016-22	24474-00	28007-05Z	35001-03	50319-30	52687-02
8023-14	8653-10GL	24060-08	24474-01	28007-05ZA	35001-04	50319-40	52687-03
8023-15	8653-10SU	24060-10	24474-02	28007-06	35001-05	50320-20	52687-04
8023-16	8653-11	24060-12	24474-03	28007-06A	35001-06	50320-30	52687-05
8023-17	8662-14	24060-14	24474-04	28007-06Z	35001-07	50323-26	52687-06
8023-18	8663-14	24060-16	24474-05	28007-06ZA	35001-08	50325-26	52689-01
8023-20	8672-10	24062-08	24474-06	28008-01	35001-09	50328-24	52689-02
8024-10	8673-10	24062-10	24480-11	28008-01A	35001-10	50336-28	52689-03
8024-11	8674-10	24062-12	24480-22	28008-01Z	35001-11	50338-28	52689-04
8024-13	8687-10	24062-14	24480-33	28008-01ZA	35001-12	50340-28	52689-05
8024-14	8688-10	24062-16	24485-00	28008-02	35001-20	50342-28	52689-06
8024-14SU	8689-10	24064-08	24485-01	28008-02A	35002-00	50345-40	52689-07
8024-15	8692-14	24064-10	24485-12	28008-02Z	35002-01	50345-50	52689-08
8024-16	8692-14SU	24064-12	24487-00	28008-02ZA	35002-02	50347-01	52689-09
8024-17	8693-14	24064-14	24487-01	28008-03	35002-03	50347-02	52689-10
8024-18	8693-14-B	24064-16	24487-12	28008-03A	35002-04	50351-01	52691-06*
8024-20	8693-14SU	24074-04*	24489-01	28008-03Z	35002-05	50424-01	52696-01
8025-10	8694-09	24077-01	24489-10	28008-03ZA	35002-06	50424-02	52734-41
8025-11	8696-11	24077-02	24489-20	28008-04	35002-07	50424-03	52742-25
8025-13	8697-10	24077-03	24490-01	28008-04A	35002-08	50424-04	52742-40
8025-14	8697-12	24077-04	24490-10	28008-04Z	35002-09	50424-10	52742-65*
8025-14L	8697-12SU	24077-05	24490-20	28008-04ZA	35002-10	50424-11	52744-23*
8025-15	8697-15	24080-06	24549-03	28008-05	35002-11	50424-12	52744-24
8025-16	8697-15SU	24080-09	24563-19	28008-05A	35002-12	50424-13	52746-19
8025-17	8697-15SUB	24080-13	24568-17	28008-05Z	35002-20	50424-14	52747-19
8025-18	8697-18SU	24080-16	24569-17	28008-05ZA	35003-00	50425-01	52770-61
8025-20	8697-20SU	24080-19	24570-17	28008-06	35003-01	50425-02	52776-60
8027-14	8697-20SUB	24080-25	24572-01	28008-06A	35003-02	50425-03	52777-40
8027-16	8697-23SUB	24080-32	24572-02	28008-06Z	35003-03	50425-04	52778-60
8027-18*	8698-12	24081-06	24575-01	28008-06ZA	35003-04	50425-52*	53361-30
8028-14	8698-15	24081-09	24575-02	28009-01A	35003-05	50426-02	53361-40
8028-16	8699-11SU	24081-13	24575-03	28009-01Z	35003-06	50426-05	53361-50
8028-18*	8702-20	24081-16	24575-04	28009-01ZA	35003-07	50426-06	53361-60
8029-14	8702-23	24081-19	24575-05	28009-02A	35003-08	50426-07	53364-30
8029-16	8703-20	24081-25	24578-40	28009-02Z	35003-09	50426-08	53364-40
8029-18*	8703-20-C	24081-32	24578-45	28009-02ZA	35003-10	50426-42	53364-50
8033-11	8703-23	24082-06*	24578-55	28009-03A	35003-11	50427-01	53364-60
8033-13	8704-18SU	24082-09	24582-23	28009-03Z	35003-12	50427-02	53365-30
8033-14	8704-23SU	24082-12*	24588-20	28009-03ZA	35003-20	50427-03	53365-40
8033-15	8705-18SU	24082-13*	24601-15	28009-04A	35004-00	50427-04	53365-50





8033-16	8705-23SU	24082-14*	24602-15	28009-04Z	35004-01	50432-01	53365-60
8033-17	8706-15SU	24082-16	24603-17	28009-04ZA	35004-02	50432-02	53367-30
8033-18	8707-15SU	24082-19	24604-17	28009-05A	35004-03	50455-08	53367-40
8033-20	8707-15SUB	24082-25	24605-10*	28009-05Z	35004-04	50457-08	53367-50
8034-11	8708-15	24083-09	24605-11*	28009-05ZA	35004-05	50459-08	53367-60
8034-13	8722-14	24083-13*	24605-12	28009-06A	35004-06	50459-14	53410-02
8034-14	8724-17	24083-16	24605-14	28009-06Z	35004-07	50459-16	53410-03
8034-14SU	8725-17	24083-19	24606-01	28009-06ZA	35004-08	50461-08	53410-04
8034-15	8727-11	24083-25	24606-02	28010-01	35004-09	50461-14	53410-05
8034-16	8728-18	24084-05	24606-03	28010-01A	35004-10	50463-09	53410-06
8034-17	8728-21	24084-10	24607-01	28010-01ZA	35004-11	50467-08	53410-08
8034-18	8728-24	24084-16	24607-02	28010-02	35004-12	50469-01	53410-09
8034-20	8730-21*	24084-25	24608-19	28010-02A	35004-13	50469-02	53410-12
8035-11	8732-19	24085-05	24609-20	28010-02ZA	35004-14	50469-03	53410-13
8035-13	8740-15	24085-10	24610-18*	28010-03	35004-15	50472-08	53410-14
8035-14	8741-15	24085-16	24610-21	28010-03A	35004-16	50472-45	53410-16
8035-15	8742-15	24085-25	24611-19	28010-04	35004-20	50472-90	53550-20
8035-16	8744-22	24086-04	24612-19	28010-04A	35005-00	50491-01	53581-01*
8035-17	8745-16	24086-06	24612-20*	28010-05	35005-01	50491-02	54112-28
8035-18	8745-23	24086-08	24613-00	28010-05A	35005-02	50491-03	54113-28
8035-20	8745-28	24086-10	24613-02	28010-06	35005-03	50493-19	54114-13
8037-14	8746-16	24086-12	24613-04*	28010-06A	35005-04	50494-14	54115-22
8037-16	8746-23	24087-04	24614-01	28011-01	35005-05	50538-00	54118-12
8038-14	8746-28	24087-06	24614-02	28011-01A	35005-06	50539-01	54420-13
8038-16	8747-16	24087-08	24614-18	28011-01Z	35005-07	50550-01	54421-13
8038-18*	8747-19	24087-10	24614-20	28011-01ZA	35005-08	50550-02	56100-04 W
8039-14	8747-23	24087-12	24615-18	28011-01Z-C*	35005-09	50550-03	56155-01 W
8039-16	8747-28	24090-04	24616-01+	28011-02	35005-10	50550-04	56155-02 W
8039-18*	8752-19	24090-06	24616-04	28011-02A	35005-11	50550-06	56158-01 W
8043-13	8753-19	24090-08	24616-05	28011-02-C	35005-12	50552-01	56158-02 W
8043-14	8754-16	24090-10	24616-23*	28011-02Z	35005-13	50552-03	56194-13 W
8044-13	8754-19	24090-12	24617-23	28011-02ZA	35005-14	50552-04	56195-02 W
8044-14	8755-16	24090-14	24618-23	28011-02Z-C*	35005-15	50552-05	56196-12 W
8045-13	8755-19	24091-04	24619-04	28011-03	35005-16	50552-07	56730-03 W
8045-14	8756-14	24091-06	24619-06	28011-03A	35005-20	50572-01	56740-03 W
8053-13	8756-16	24091-08	24619-08	28011-03-C	35006-00	50572-02	607-8112
8053-14	8756-19	24091-10	24619-23+	28011-03Z	35006-01	50573-01	607-8116
8054-13	8757-16	24091-12	24620-20	28011-03ZA	35006-02	50573-02	610-8100
8054-14	8757-19	24091-14	24621-02	28011-03Z-C*	35006-03	50574-01	610-8102
8055-13	8758-16	24092-04*	24621-03*	28011-04	35006-04	50574-02	610-8106
8055-14	8758-19	24092-06	24621-04	28011-04A	35006-05	50575-01	610-8110
8056-13*	8762-24	24092-08	24621-06	28011-04-C	35006-06	50575-02	610-8114
8062-14	8767-23	24092-10	24621-08	28011-04Z	35006-07	50579-01	610-8200
8063-14	8769-22	24092-12*	24621-13	28011-04ZA	35006-08	50579-02	610-8201





8072-14	8771-12	24092-13	24621-15	28011-04Z-C*	35006-09	50579-03	610-8202
8073-14	8771-25	24092-15	24621-18	28011-05	35006-10	50579-04	610-8206
8082-14	8771-45	24092-18*	24622-04	28011-05A	35006-11	50579-05	703-205LL
8102-15	8771-45SU*	24092-20	24622-18	28011-05-C	35006-12	50610-01	704-140
8103-15	8771-60	24092-25	24624-22	28011-05Z	35006-13	50610-02	704-145
8103-15SU	8772-12	24092-30	24625-22	28011-05ZA	35006-14	50610-10	704-146
8112-15	8772-25	24092-38	24627-17	28011-05Z-C*	35006-15	50610-11	708-400LL*
8113-15	8772-45	24092-44	24628-36*	28011-06	35006-16	50610-12	709-103*
8122-16	8772-60	24092-50	24629-23	28011-06A	35006-20	50610-13	712-442
8122-16SU	8772-90	24093-04*	24630-23	28011-06Z	35007-00	50610-14	712-452
8122-21	8773-25	24093-06	24630-27*	28011-06ZA	35007-01	50610-15	712-462
8123-16	8773-45	24093-08	24631-16	28012-01	35007-02	50611-01	712-472
8123-16-B	8773-60	24093-10	24632-02	28012-02	35007-03	50611-02	779-1001
8123-16-D	8774-12	24093-13	24632-03	28012-03	35007-04	50611-10	783-720
8123-16SU	8774-22*	24093-15	24632-04*	28012-04	35007-05	50611-11	783-721
8123-21	8774-25	24093-18*	24633-01	28012-05	35007-06	50611-12	790-100
8125-16-C	8774-45	24093-20	24633-02	28012-06	35007-07	50611-13	790-102
8125-16-D	8774-60	24093-25	24633-03*	28013-01	35007-08	50611-14	790-104
8125-16-E	8774-90	24093-30	24633-16*	28013-01A	35007-09	50611-15	790-106
8132-13	8775-23	24093-38	24634-24	28013-01Z	35007-10	50614-01	790-120
8132-14	8776-23	24093-44	24635-14	28013-02	35007-11	50614-02	790-196
8134-11	8777-16	24093-50	24635-14D	28013-02A	35007-12	50614-10	790-197
8134-11SU	8777-23	24094-02*	24636-14*	28013-02Z	35007-13	50614-11	790-198
8135-11	8832-18	24094-03	24637-17	28013-02ZA	35007-14	50614-12	790-200
8135-11SU	8837-17	24094-04	24639-17	28013-02Z-C	35007-15	50614-13	790-201
8138-16	8842-18	24094-06	24641-15	28013-03	35007-16	50614-14	790-202
8138-21	8843-18	24094-08	24642-16*	28013-03A	35007-20	50614-15	790-203
8139-16	8847-17	24094-10	24643-16	28013-03Z	35012-27*	50615-01	790-204
8139-16SR	8849-23	24094-12	24645-16	28013-03ZA	35015-24*	50615-02	790-205
8139-21	8852-17	24094-14*	24646-23*	28013-03Z-C	35016-01	50615-10	790-216
8142-14	8857-17	24094-15	24646-27	28013-04	35016-03	50615-11	790-216-34*
8152-14	8857-17SU	24094-34	24647-27	28013-04A	35510-01	50615-12	790-217
8152-16	8859-17	24094-36	24649-19	28013-04Z	35510-02	50615-13	790-217-34
8152-18	8862-19	24095-03	24650-14	28013-05	35510-03	50615-14	790-218
8154-14	8867-19	24095-04	24650-17	28013-05A	35510-04	50615-15	790-218-34*
8154-18	8872-19	24095-06	24650-19	28013-05Z	35510-05	50620-13	790-220
8155-14	8882-13*	24095-08	24650-22	28013-05ZA	35510-06	50621-01	790-221
8155-18	8882-30	24095-10	24651-14	28013-06	35511-01	50621-04	790-222
8188-11	8912-27	24095-12	24651-17	28013-06A	35511-02	50621-05	790-223
8188-11SU	8922-21	24095-15	24651-19	28013-06Z	35511-03	50621-20	790-224
8188-14	8927-21	24096-07	24651-22	28014-02*	35511-04	50621-21*	790-225
8188-14SU	8932-21	24096-10	24653-15*	28014-05*	35511-05	50621-22	790-236
8189-11	8935-30*	24096-15	24654-15	28015-01A	35511-06	50622-16*	790-237
8189-11SU	8936-30*	24096-20	24655-15	28015-01Z	35512-01	50622-17	790-238





8189-11SUB	8939-25	24096-25	24660-18	28015-02A	35512-02	50623-12	790-240
8189-14	8940-31*	24096-30	24661-18	28015-02Z	35512-03	50624-12	790-241
8189-14SU	9022-13	24097-07	24666-24*	28015-02ZA	35512-04	50625-01*	790-242
8193-14	9022-16	24097-10	24667-24*	28015-03A	35512-05	50625-02*	790-243
8197-14 SU TC	9032-16	24097-15	24670-27	28015-03Z	35512-06	50625-04*	790-244
8198-14	9042-10	24097-20	24671-23*	28015-04A	35513-01	50627-18*	790-245
8198-14SU	9072-14	24097-25	24671-27	28015-04Z	35513-02	50628-16	790-251
8198-14SUB	9077-20*	24097-30	24673-22*	28015-05A	35513-03	50630-00	790-255
8250-10	9077-22	24098-05*	24673-27	28015-05Z	35513-04	50630-01	790-256
8250-14*	9081-23	24098-10*	24675-27	28015-06A	35513-05	50630-02	790-258
8272-14	9083-22	24098-15*	24677-24	28015-06Z	35513-06	50630-10	790-260
8272-14SU	9085-24*	24099-05*	24680-24	28016-03	35730-01	50630-20	790-262
8272-14SUB	9085-25*	24099-10*	24681-32	28016-05	35730-02	50634-01	790-263
8272-15	9091-23	24099-15*	24682-32	28017-03	35730-03	50634-02	790-264
8272-15SU	9093-22	24100-08	24683-36	28017-05	35730-13*	50634-03	790-368
8272-17	9102-15*	24100-10	24685-24	28018-03	35732-30	50690-13	790-369
8272-17L*	9107-18	24100-12	24685-26	28018-05	35735-02	50705-17	790-370
8272-17SU	9132-17	24100-15	24687-26	28019-03	38110-35	50710-15	790-385
8272-17SUB	9134-11	24100-20	24688-27	28019-05	38112-08	50712-16	790-387
8272-23	9134-11 SU TC	24100-25	24693-23	28022-03	38112-10	50714-16	790-410
8272-23SU	9135-11	24102-08	24694-34	28022-05	39100-01	50753-16	790-412
8272-28*	9135-11 SU TC	24102-10	24851-01	28023-03	39100-02	50755-16	790-667
8272-40*	9136-13*	24102-12	24851-02	28023-05	39100-04	50757-16	790-668
8273-14	9137-18	24102-15	24851-03	28024-01	39100-05	50761-18	790-669
8273-14SU	9138-20*	24102-18*	24851-04	28024-02	39100-06	50763-18	790-670
8273-14SUB	9142-19	24102-20	24854-05	28025-02	39101-01	50766-04	790-671
8273-15	9152-13	24102-25	24854-09	28025-03	39101-02	50766-05	790-672
8273-15SU	9152-13-B	24103-06	24855-01	28025-04	39101-03	50766-06	790-673
8273-15SUB	9152-13SU	24103-08	24855-02	28025-05	39101-04	51120-01	790-674
8273-16SU	9153-13*	24103-10	24855-03	28026-02	39101-05	51120-02	790-681
8273-17	9153-13-C	24103-12	24855-04	28026-03	39178-01	51121-01	790-682
8273-17L*	9157-16	24103-15	24855-05	28026-04	39178-02	51121-02	790-683
8273-17SU	9162-19	24103-20	24855-06	28027-02	39201-01	51123-01	791-092*
8273-17SUB	9163-20SU	24103-25	24855-07	28027-03	39202-10	51123-02	791-093*
8273-18L	9163-20SUB	24106-08	24855-08	28027-04	39202-11	51123-03	791-097
8273-23	9164-20	24106-10	24855-09	28028-03*	39205-38	51123-04	791-3514
8273-23SU	9165-20 SU TC	24106-12	24855-10	28028-13*	39400-00	51124-01	791-3515
8273-28*	9167-19	24106-15	24855-11	28030-13	39400-01	51124-02	791-3516
8273-40*	9168-19	24106-20	24855-12	28030-13Z	39400-02	51124-03	791-3517
8275-22	9168-19SU	24106-25	24864-14	28030-15	39400-03	51131-12	793-180
8276-22	9169-19-B	24107-08	24898-14	28030-18	39400-04	51133-12	793-181
8278-14	9169-19SU	24107-10	24899-21	28030-18Z	39400-05	51141-10	793-182
8278-17	9170-19	24107-12	24907-15	28030-20Z-C	39400-06	51143-12	793-183
8278-17L*	9170-23	24107-15	24908-20	28030-22	39400-08	51144-12*	793-184





8278-23	9171-19	24107-20	24909-21	28031-13	39400-09	51144-13*	793-185
8279-14	9171-19-C	24107-25	24915-30	28031-13Z	39400-10	51147-12	793-364
8279-17	9171-23	24111-02	24915-40	28031-15	39400-11	51155-01	793-365
8279-17SR	9172-18	24111-03	24915-50	28031-18	39400-14	51155-02	793-366
8279-23	9172-19 SU TC	24111-04	24915-60	28031-18Z	39400-15	51155-03	793-367
8281-17*	9173-18	24111-05	24915-70	28031-22	39400-16	51155-04	793-368
8282-15	9173-19 SU TC	24111-06	24916-00	28033-13	39400-18	51155-05	793-369
8282-15SU	9174-14	24111-07	24917-30	28033-13Z	39400-20	51156-01	793-376
8282-17	9174-14SU	24111-08	24917-40	28033-15	39400-21	51156-02	793-377
8282-17-E*	9174-18	24113-02	24917-50	28033-18	39400-22	51156-03	794-385
8282-17L*	9174-18SU	24113-03	25303-13	28033-18Z	39400-24	51156-04	794-387
8282-17SU	9175-14	24113-04	25303-19	28033-22	39400-26	51156-05	794-390
8282-19	9175-14-C	24113-05	25304-20	28034-13	39400-28	51156-06	794-391
8282-21	9175-14SU	24113-06	25360-24	28034-13Z	39400-30	51156-21	794-394
8283-15	9175-18	24113-07	25372-22	28034-15	39400-31	51156-22	794-395
8283-15SU	9175-18SU	24113-08	25379-20	28034-18	39400-32	51156-23	AH 100-04
8283-16+	9176-23	24117-02	25380-29	28034-18Z	39400-33	51156-24	AH 110-03
8283-16SU	9176-23SU	24117-03	25386-21	28034-20	39400-34	51156-25	AJ 023-11
8283-17	9176-23SUB	24117-04	25387-21	28034-20Z-C*	39400-35	51156-26	AJ 023-13
8283-17-A	9177-23	24117-05	25390-25	28034-22	39400-40	51157-01+	AJ 023-14
8283-17-B	9177-23SU	24117-06	25395-01	28035-13	39400-43	51157-02+	AJ 023-15
8283-17-D	9177-23SUB	24117-07	25395-02	28035-13Z	39400-48	51157-03+	AJ 023-16
8283-17-E	9178-14 SU TC	24117-08	25395-03	28035-15	39400-50	51157-04+	AJ 023-17
8283-17SU	9178-18 SU TC	24211-02	25395-04	28035-18	39400-51	51157-05+	AJ 023-18
8283-19	9178-23 SU TC	24211-03	25395-05	28035-18Z	39400-52	51158-01	AJ 024-11
8283-19SU	9179-11	24211-04	25395-06*	28035-20	39400-53	51158-02	AJ 024-13
8283-21	9179-14 SU TC	24211-05	25410-15	28035-22	39400-54	51158-03	AJ 024-14
8283-23SU	9179-18 SU TC	24211-06	25410-19	28037-13	39400-55	51158-04	AJ 024-15
8288-15	9179-23 SU TC	24211-07	25411-16	28037-13Z	39400-56	51159-01	AJ 024-16
8288-17	9180-11	24211-08	25411-17	28037-15	39400-57	51159-02	AJ 024-17
8289-15	9181-14SU	24213-02	25411-18	28037-18	39400-60	51159-03	AJ 024-18
8289-17	9182-17	24213-03	25411-19	28037-18Z	39400-61	51160-01	AJ 025-11
8292-14	9182-17SU	24213-04	25435-01	28037-20	39400-67	51160-02	AJ 025-13
8292-17	9182-17SUB	24213-05	25435-02	28037-22	39400-75	51160-03	AJ 025-14
8293-14	9182-18SU	24213-06	25438-04*	28038-13	39400-76	51160-04	AJ 025-15
8293-17	9182-18SU TC	24213-07	25439-01*	28038-13Z	39400-78	51160-05	AJ 025-16
8295-14	9182-20	24213-08	25439-02*	28038-15	39400-79	51161-01	AJ 025-17
8295-22	9182-20SU	24217-02	25500-18	28038-18	39400-80	51161-02	AJ 025-18
8323-11	9182-20SUB	24217-03	25501-18	28038-18Z	39400-81	51161-03	AJ 033-11
8323-11SU	9183-17	24217-04	25502-18	28038-22	39400-85	51161-04	AJ 033-13
8323-14	9183-17SU	24217-05	25520-16	28039-13	39400-86	51162-01	AJ 033-14
8323-14-E	9183-17SUB	24217-06	25520-18	28039-13Z	39400-90	51162-03	AJ 033-15
8323-14SU	9183-18 SU TC	24217-07	25521-16	28039-15	39400-91	51165-01	AJ 033-16
8324-14	9183-18SU	24217-08	25521-18	28039-18	39400-93	51165-02	AJ 033-17





8325-14	9183-20	24220-03	25522-16	28039-18Z	39400-94*	51165-03	AJ 033-18
8327-11	9183-20SU	24220-05	25522-18	28039-22	39400-95	51165-04	AJ 034-11
8327-14	9183-20SUB	24220-06	25523-18	28041-13	39400-97	51166-01	AJ 034-13
8328-14	9184-16SU	24220-08	25524-18	28041-13Z	39400-98	51171-11	AJ 034-14
8329-14	9185-23SU	24220-09	25526-16	28041-18	39401-02	51171-12	AJ 034-15
8333-11	9186-23SU	24221-06	25526-18	28041-18Z	39401-03	51171-13	AJ 034-16
8333-11-A	9187-12SU	24221-09*	25527-16	28041-22	39401-04	51173-12	AJ 034-17
8333-11-B	9188-11	24221-13	25527-18	28042-18	39401-05	51173-21*	AJ 034-18
8333-11-D	9188-11 SU TC	24221-19	25550-19	28044-21*	39401-68*	51173-31	AJ 035-11
8333-11-E	9188-14	24221-25	25550-22	28050-02	39401-69*	51174-01	AJ 035-13
8333-12SU*	9188-14 SU TC	24221-32	25551-22	28050-03	39401-70	51174-02	AJ 035-14
8333-14	9189-11	24222-06	25556-18	28050-04	39401-71	51175-01	AJ 035-15
8333-14-A	9189-11 SU TC	24222-13	25556-22	28050-06	39401-72	51175-02	AJ 035-16
8333-14-B	9189-14	24222-19	25557-18	28051-02	39401-73	51175-03	AJ 035-17
8333-14-E	9189-14 SU TC	24222-25	25557-22	28051-03	39402-10	51175-05	AJ 035-18
8333-14SU	9190-10	24222-32	25558-22	28051-04	39402-12	51176-01	AK 432-11
8333-14SUB	9191-18	24223-06	25559-22	28053-02	39402-20	51176-11	AK 434-11
8334-14	9191-18SU	24223-13	25600-15	28053-03	39402-25	51184-20	AK 436-11
8335-14	9191-18SUB	24223-19	25601-15	28053-04	39402-30	51186-21*	AK 438-11
8337-11	9192-17	24223-25	25602-01	28070-19	39402-32	51194-13	Z0000025127
8337-14	9192-20	24223-32	25602-02	28074-25	39402-34	51195-13	Z0000087375
8337-14SR	9193-14	24224-06	25602-03	28075-25	39403-10	51196-12	Z0000093030
8338-14	9193-17	24224-13	25602-04	28080-02	39404-10	51197-12	Z0000093216
8339-14	9193-20	24224-19	25602-18	28080-03	39599-01	51198-04	Z0000093552
8343-15	9194-18	24224-25	25603-01	28080-04	39599-02	51230-02	Z0000094794
8343-18	9195-14	24224-32	25603-02	28081-02	39599-03	51230-04	Z0000095290
8343-18-B	9196-14	24225-06	25603-03	28081-03	39599-04	51230-05	Z0000097731
8343-18SU	9197-14 SU TC	24225-13	25603-04	28081-04	39599-05	51230-06	Z0000100289
8343-18SUB	9198-11	24225-19	25603-18	28081-05	39599-06	51245-20	Z0000105016
8343-20	9198-11SU	24225-25	25606-18	28081-06	39599-08	51250-20	Z0000105079
8343-20-E	9198-14	24225-32	25607-18	28082-02	39599-09	51251-00	Z0000106100
8343-20SU	9198-14SU	24227-04	25610-01	28082-03	39599-10	51251-01	Z0000106718
8343-23	9199-11	24227-06	25610-02	28082-04	39599-11	51251-02	Z0000109527
8343-23SU	9199-11SU	24228-03	25610-03	28082-05	39599-12	51251-03	Z0000109678
8343-25	9199-13	24228-04	25610-04	28082-06	39599-16	51302-01	Z0000109680
8343-28	9200-11	24228-05	25610-15	28083-03	39599-17	51302-02	Z0000109682
8343-30	9200-11 SU TC	24228-06	25611-01	28083-04	39599-18	51303-01	Z0000110039
8344-18	9201-11	24239-11	25611-02	28084-03	39599-19S	51303-02	Z0000110416
8344-20	9201-11 SU TC	24239-12	25611-03	28084-04	39599-19W	51304-01	Z0000113442
8344-23	9202-13	24248-03	25611-04	28084-05	39600-00	51304-02	Z0000114383
8344-25	9203-10	24248-05	25611-15	28085-03	39600-01	51305-01	Z0000114410
8344-28	9204-11	24248-08	25612-12	28085-04	39600-02	51306-15	Z0000115910
8344-30	9205-18	24249-12	25613-12	28085-05	39600-03	51308-15	Z0000115925
8345-18	9206-11	24249-14	25616-12	28086-03	39600-04	51313-13	Z0000115934





8345-20	9207-11	24249-16	25617-12	28086-04	39600-05	51315-23	Z0000115987
8345-23	9208-12	24250-02	25620-18	28087-03	39600-06	51317-22	Z0000115994
8345-25	9208-12SU	24250-03	25621-18	28087-04	39600-08	51321-21+	Z0000115999
8345-28	9208-12SUB	24250-04	25670-14*	28088-03	39600-09	51323-23	Z0000116003
8345-30	9208-16	24250-05	25670-15	28088-04	39600-10	51330-11	Z0000117905
8347-18	9209-12SU	24250-06	25671-14*	28089-01	39600-11	51330-12	Z0000117997
8347-20	9209-16	24250-08	25671-15	28089-02	39600-12	51335-03	Z0000118030
8347-23	9210-16	24250-10	25672-15	28089-03	39600-16	51342-12	Z0000118203
8347-25	9211-09	24251-15	25673-15	28090-01	39600-17	51344-15	Z0000118633
8347-28	9211-15*	24252-16	25675-15	28090-01Z	39600-18	51346-15	Z0000118634
8347-30	9211-15SU*	24253-04	25677-18	28090-02	39600-19S	51350-01	Z0000119422
8348-18	9211-18SU*	24253-06	25678-15	28090-02Z	39600-19W	51355-01	Z0000119464
8348-20	9212-09	24253-08	25679-15	28090-02Z-C*	39600-20	51361-02	Z0000119590
8348-23	9212-13*	24254-03	25680-15	28090-03	39600-22	51364-21	Z0000119710
8349-18	9213-09	24255-04	25681-15	28090-03Z	39600-23	51366-01	Z0000119714
8349-20	9214-14SU	24255-06	25715-15	28090-04	39600-24	51366-02	Z0000119742
8349-23	9215-10SU	24255-08	25750-01	28090-04Z	39600-25	51366-03	Z0000119746
8351-20*	9216-10SU	24261-16	25750-02	28090-05	39600-26	51366-04	Z0000119844
8351-23*	9217-16SU	24263-16	25750-03	28090-05Z	39600-27	51366-05	Z0000119923
8351-28*	9222-01	24270-07	25751-01	28090-06	39600-28	51366-06	Z0000120057
8352-18	9222-02	24270-12	25751-02	28090-06Z	39600-29	51368-12	Z0000120058
8352-18SU	9222-03	24272-19	25751-03	28091-01Z	39600-30	51369-00*	Z0000120229
8352-18SUB	9223-01	24277-02	25751-04	28091-02Z	39650-01	51405-01	Z0000120458
8353-15	9223-02	24277-06	25751-05	28091-03Z	39650-02	51405-02	Z0000120475
8353-15-A	9223-03	24278-04	25751-06	28091-04Z	39650-20	51405-03	Z0000120533
8353-15-B	9230-12	24278-05	25752-01	28091-05Z	39651-01	51405-04	Z0000120534
8353-16	9243-11	24278-06	25752-02	28091-06Z	39651-02	51405-05	Z0000120553
8353-18	9244-11	24278-07	25752-03	28094-01	39651-03	51410-03	Z0000120561
8353-18-A	9245-11	24278-08	25752-04	28094-02	39651-04	51410-04	Z0000120595
8353-18-B	9252-05	24279-10	25753-01	28094-03	39651-05	51410-05	Z0000120717
8353-18-C	9257-05	24279-12	25753-02	28094-04	39651-06	51411-03	Z0000121458
8353-18-E	9262-09*	24279-14	25753-03	28094-05	39651-07	51411-04	Z0000121460
8353-18SU	9262-11	24279-16	25753-04	28094-06	39651-08	51411-05	Z0000121623
8353-18SUB	9263-11	24280-04	25753-05	28095-01Z	39651-09	51412-01	Z0000121794
8353-20	9272-09	24280-06*	25753-06	28095-02	39651-10	51412-02	Z0000122120
8353-20-B	9282-08	24280-07	25754-12	28095-02Z	39653-02	51412-03	Z0000122140
8353-20-C	9282-08S*	24280-08*	25754-25	28095-02Z-C	39653-08	51412-04	Z0000122141
8353-20-E	9283-08	24280-09	25754-45	28095-03	39653-10	51412-05	Z0000122232
8353-20SU	9300-14	24280-12	25754-60	28095-03Z	39654-01	51412-06	Z0000122346
8353-23	9319-14*	24281-04	25754-90	28095-03Z-C*	39654-02	51413-01	Z0000122348
8353-23-B	9322-10	24281-07	25755-12	28095-04	39654-03	51413-02	Z0000122684
8353-23-E	9332-10	24281-09	25755-25	28095-04Z	39654-04	51413-03	Z0000122691
8353-23LH*	9342-11	24281-12	25755-45	28095-04Z-C*	39654-05	51413-04	Z0000122838
8353-23SU	9347-11	24283-06	25755-60	28095-05	39654-06	51413-05	Z0000123004





8353-23SUB	9352-11	24284-16	25755-90	28095-05Z	39654-07	51415-01	Z0000123218
8353-25	9357-11	24286-16	25756-12	28095-06	39654-08	51415-02	Z0000123253
8353-25-E	9358-14*	24287-03	25756-25	28095-06Z	39654-09	51415-03	Z0000123370
8353-25SU	9362-11	24287-04	25756-45	28095-06Z-C*	39654-10	51415-04	Z0000123425
8353-28	9362-22	24287-07	25756-60	29092-23	39654-11	51513-01	Z0000123708
8353-30	9367-11	24287-09	25756-90	29092-33	39654-12	51513-02	Z0000123818
8354-18	9367-22	24287-12	25757-12	29092-40	39654-15	51513-03	Z0000123842
8354-20	9422-12	24290-10	25757-25	29094-23	39654-16	51522-01	Z0000124042
8354-23	9423-12	24290-12	25757-45	29094-33	39657-04	51522-02	Z0000124098
8354-25	9432-11	24290-14	25757-60	29094-40	39658-02	51522-03	Z0000124680
8354-28	9432-11-B	24292-03	25757-90	29098-23	39658-04	51655-20	Z0000124814
8354-30	9432-11-E	24292-04	25758-12	29098-33	39659-02	51657-20	Z0000125120
8356-18	9432-11GL	24294-03*	25758-25	29098-40	39659-04	51661-01	Z0000125793
8356-20	9432-11SU	24294-04*	25758-30*	29100-28	39660-02	51661-02	Z0000125811
8356-23	9432-11SU*	24294-07*	25758-31*	29100-35	39661-02	51661-03	Z0000125966
8356-25	9432-11SUB*	24294-12*	25758-32*	29100-42	39670-01	51663-18	Z0000126703
8356-28	9434-11	24295-02	25758-33*	29101-28	39670-02	51664-16	Z0000126736
8356-30	9434-11-B	24295-03*	25758-34*	29101-35	39670-03	51664-17	Z0000126866
8357-18	9434-11-C	24295-04	25758-45	29101-42	39670-04	51665-21	Z0000126925
8357-18SR	9434-11GL	24295-05*	25758-60	29102-28	39671-02	51667-19	Z0000126926
8357-18SUTC*	9434-11SR	24295-06	25758-90	29102-35	39671-04	51667-20	Z0000126927
8357-18WB	9434-11SU	24295-08	25759-12	29102-42	39673-02	51669-21	Z0000127196
8357-20	9436-11	24295-10	25759-25	29103-28	39673-04	51670-01	Z0000127214
8357-20SR	9436-11D	24295-12	25759-45	29103-35	39674-04	51670-02	Z0000127678
8357-20W*	9438-11	24295-15	25759-60	29103-42	39680-01	51671-21	Z0000127706
8357-23	9452-10	24295-18	25759-90	29129-35	39680-04	51672-21	Z0000127798
8357-23L	9452-10SU	24296-01	25760-01	30258-34	39680-05	51673-18*	Z0000128030
8357-23SR	9452-11	24296-02	25760-02	30285-01	39680-06	51675-18	Z0000128094
8357-25	9452-11SU	24296-03	25760-03	30285-02	39681-01	51677-18	Z0000128269
8357-25SR	9452-12*	24298-01	25760-04	30285-03	39681-02	51683-01	Z0000128435
8357-28	9452-12SU*	24298-02	25760-05	30285-04	39681-03	51683-02	Z0000128673
8357-30	9453-10	24298-03	25760-06	30285-05	39681-04	51683-03	Z0000128674
8357-35	9453-10SU	24390-01	25760-07	30285-06	39681-05	51685-21	Z0000128780
8358-18	9453-11	24390-02	25760-08	30285-07	39681-06	51686-01	Z0000128805
8358-20	9453-11-B	24390-03	25761-01	30285-08	39682-01	51686-03	Z0000128823
8358-23	9453-11SU	24390-04	25761-02	30290-28	39682-02	51686-05	Z0000128898
8359-18	9453-11SUB	24390-05	25761-03	30291-28*	39682-03	51686-07	Z0000128906
8359-20	9453-12SU*	24390-06	25761-04	30295-21	39682-04	51686-09	Z0000129021
8359-23	9454-10	24390-10	25761-05	30295-22	39682-05	51700-01	Z0000129248
8361-16*	9455-10	24390-20	25761-06	30295-23	39682-06	51700-02	Z0000129288
8363-14	9456-10	24390-30	25761-07	30423-00	39689-02	51700-03	Z0000129351
8365-14	9457-10	24391-01	25761-08	30424-00	40000-01	51710-03	Z0000129352
8366-14	9462-10	24391-02	25769-45	30424-01	40000-03	51710-04	Z0000129353
8367-14	9462-10SU	24391-03	25856-22	30424-02	40000-04	51711-03	Z0000129354





8369-16W	9462-11	24391-04	25857-22	30424-03	40002-02	51711-04	Z0000129379
8370-16*	9462-11-B	24391-05	25902-18+	31100-02	40004-01	51712-03	Z0000129535
8370-18	9462-11-B	24391-06	25903-15*	31100-12	40006-02	51715-02	Z0000129536
8370-20	9462-11SU	24391-10	28001-01	31100-16	40006-05	51716-02	Z0000129769
8370-23	9463-10	24391-20	28001-01A	31100-18	40007-01	51720-03	Z0000129809
8370-28*	9463-10SU	24391-30	28001-01ZA	31100-23	40008-02	51722-03	Z0000129861
8371-16*	9463-10SUB	24392-01	28001-02	31101-12	40009-03	51724-03	Z0000130016
8371-18	9463-11	24392-10	28001-02A	31101-16	40010-01	51725-04	Z0000130152
8371-18SR	9463-11-B	24392-20	28001-02ZA	31101-18	40010-02	51727-03	Z0000130159
8371-20	9463-11SU	24392-30	28001-03	31101-23	40011-01	51728-04	Z0000130240
8371-23	9463-11SUB	24393-01	28001-03A	31102-12	40011-02	51729-00	Z0000130291
8371-25*	9463-13SU*	24393-10	28001-03ZA	31102-16	40011-03	51729-01	Z0000130323
8371-28*	9469-10	24393-30	28001-04	31102-18	40011-04	51729-02	Z0000130369
8372-14	9470-10*	24394-10	28001-04A	31102-23	40011-05	51729-03	Z0000130370
8372-14SU	9471-09	24394-20	28001-04ZA	31103-12	40011-06	51729-04	Z0000130393
8373-14	9472-10	24394-30	28001-05	31103-16	40011-07	51729-05	Z0000130402
8373-14-A	9473-10	24394-40	28001-05A	31103-18	40011-08	51729-06	Z0000130426
8373-14-B	9477-09	24395-10	28001-05ZA	31103-23	40011-09	51730-03	Z0000130428
8373-14-E	9482-11	24395-20	28001-06	31110-12	40011-10	51730-04	Z0000130429
8373-14SU	9482-11SU	24395-30	28001-06A	31110-16	40011-11	51731-03	Z0000130525
8373-16*	9482-11SUB	24395-40	28001-06ZA	31110-18	40011-12	51732-03	Z0000130588
8373-23*	9483-11	24395-50	28002-01	31110-23	40012-01	51733-04	Z0000130598
8374-14	9483-11SU	24396-01	28002-01A	31111-12	40012-02	51735-00	Z0000130600
8375-14	9483-11SUB	24400-13	28002-01Z	31111-16	40012-03	51735-01	Z0000130637
8380-18	9484-11	24400-14	28002-01ZA	31111-18	40013-01	51735-02	Z0000130731
8380-18SU	9485-11	24400-17	28002-02	31111-23	40013-02	51736-00	Z0000130747
8380-20	9492-11	24410-13	28002-02A	31112-12	40013-03	51737-00	Z0000130750
8380-20SU	9502-11	24410-16	28002-02-C	31112-16	40014-00	51737-01	Z0000130751
8380-23	9502-13	24410-20	28002-02Z	31112-18	40015-00	51737-02	Z0000130762
8380-23SU	9502-14	24412-21	28002-02ZA	31112-23	40016-01	51737-03	Z0000130763
8381-18	9503-11	24440-10	28002-03	31113-12	40016-02	51737-10*	Z0000130764
8381-18-B	9503-13	24440-20	28002-03A	31113-16	40018-11	51737-12*	Z0000130776
8381-18-D	9503-14	24440-30	28002-03-C	31113-18	40019-13	51740-03	Z0000130936
8381-18-E	9517-14	24440-40	28002-03Z	31113-23	40020-01	51740-04*	Z0000130939
8381-18SU	9517-17	24442-01	28002-03ZA	31120-16	40020-02	52131-13+	Z0000130947
8381-18SUB	9518-22	24442-02	28002-04	31120-23	40020-03	52601-23	Z0000130980
8381-20	9519-22	24442-03	28002-04A	31121-16	40021-01	52603-23	Z0000130981
8381-20-E	9521-17	24442-04	28002-04-C	31121-23	40023-01	52605-22	Z0000130991
8381-20SU	9522-12*	24442-05	28002-04Z	31130-16	40023-02	52607-25	Z0000131011
8381-23	9522-14	24442-06	28002-04ZA	31130-23	40023-03	52608-24*	Z0000131012
8381-23SU	9522-16*	24442-10	28002-05	31131-16	40024-00	52624-63	Z0000131062
8381-23SUB	9530-09	24442-20	28002-05A	31131-23	40027-00	52665-13	Z0000131079
8385-18	9532-09	24442-30	28002-05-C	31132-16	40029-00	52665-15	Z0000131088
8385-20	9532-10	24442-40	28002-05Z	31132-23	40040-01	52665-17	Z0000131096





8385-23	9532-11*	24443-00	28002-05ZA	31133-16	40040-02	52671-00	Z0000131114
8392-14	9532-12*	24443-01	28002-06	31133-23	40041-01	52671-01	Z0000131120
8392-14SU	9532-13	24443-02	28002-06A	31140-12*	40041-02	52671-02	Z0000131123
8394-19	9533-11	24443-03	28002-06Z	31140-16	40043-01	52671-03	Z0000131220
8423-20	9535-14	24443-04	28002-06ZA	31140-23	40043-02	52671-10	Z0000131221
8423-23	9535-18	24445-01	28003-01	31148-12*	40060-00	52673-01	Z0000131222
8424-20	9542-09	24445-02	28003-01A	31150-16	40061-00	52673-02	Z0000131223
8424-23	9542-12	24445-03	28003-01Z	31150-23	40180-02	52673-03	Z0000131303
8425-23	9552-15	24445-04	28003-01ZA	34396-23	40181-02	52673-04	Z0000131304
8433-20	9562-10	24445-05	28003-02	34397-27	40182-02	52673-05	Z0000131315
8433-23	9567-11	24445-06	28003-02A	34400-21	40185-01	52673-06	Z0000131316
8434-20	9577-10	24448-01	28003-02Z	34400-26+	40186-01	52673-07	Z0000131317
8434-23	9581-11	24448-02	28003-02ZA	34400-29	40187-01	52675-01	Z0000131318
8442-27	9593-10	24448-03	28003-03	34402-23	40188-01	52675-02	Z0000131319
8443-16*	9593-12	24448-04	28003-03A	34410-24	40189-07	52675-03	Z0000131320
8443-18*	9602-10	24448-05	28003-03Z	34420-22	40200-05	52675-04	Z0000131321
8443-27	9602-12	24448-06	28003-03ZA	34420-24	40200-10	52675-05	Z0000131322
8452-17	9603-10	24450-01	28003-04	34440-26	40200-15	52675-06	Z0000131329
8453-17	9603-12	24450-02	28003-04A	34450-28	40200-18	52675-07	Z0000131517
8462-16	9612-10	24450-03	28003-04Z	34460-23	40200-20	52675-08	Z0000135022
8462-16SR	9612-12	24450-04	28003-04ZA	34461-23	40200-25	52675-09	Z0000135030
8462-16SU	9613-10	24450-05	28003-05	34470-00	40200-30*	52675-10	Z0000138476
8462-18	9613-12	24450-06	28003-05A	34470-01	40200-35	52675-11	Z0000138486
8462-18SU	9642-12*	24450-10	28003-05Z	34470-24*	40200-40	52675-12	Z0000138490

